

NDA

MEDICAL INDUSTRIES INC.

**ELITE FAMILY OF PRODUCTS FOR
THE LONG TERM CARE MARKET**

You can trust **NDAH** for all your Acute Care needs.



Table of Contents

NOA LTC Beds

Elite Riser	4
Twin Elite Riser	6
Elite Bariatric	8
Elite R600	10
NOA Light	12

NOA Acute Care Beds

Platinum SC	20
Platinum NS	23
Platinum ML	24

Accessories

Standard Head and Foot Boards	14
Side Rails	14
Sub Acute Handset	14
Safety Lock	15
Assist Bar	15
Floor Mats	15
In Rail Pendant	16
In Bed Scale	16
Footboard Control	16
Under Bed Lighting	17
Mattresses	17
Over Bed Table	18
Deck Expansion System	18
Trapeze	19

NOA Elite Riser™

Four Year Warranty on Electronics
Lifetime warranty on bed deck and frame.

600 lb. safe working load

Faster & Quieter
The new Elite Riser swiftly moves from 7.5" to 28".

Lightweight
The Elite Riser is lighter than most high/low beds.

Self Leveling
A built-in sensor in the control box communicates with the actuators to keep the bed level at all times.

Very Mobile
The Elite Riser can be easily moved at any height and in any position.

New Safety Lock Option
Makes it one of the safest beds on the market.



Features Include:

- Bed goes straight up and straight down
- Swivel braking casters at all four corners
- Bed rolls and locks at any height
- Entrapment compliant with properly sized mattress
- Easy to maintain
- Faster
- Quieter

Optional Features:

- Safety lock for immobilization
- Head and Foot Boards - Light Oak, Dark Oak, Wild Cherry, Walnut, Mahogany
- Soft Touch Side Rails
- Sub-acute handset
- Universal Assist Bar
- Floor Pad, 70" x 23.5" x 1"
- Therapeutic Foam Mattresses
- Footboard Mounted Staff Control
- Full Overbed Trapeze
- Deck Expansion System

ITEM	76" length	80" length	84" length
bed width	35"	35"	35"
bed width, with side rails	40.0"	40.0"	40.0"
bed width, with assist bar	36.5"	36.5"	36.5"
overall bed length with head and foot boards	80.5"	84.5"	88.5"
deck height range	7.5" to 28"	7.5" to 28"	7.5" to 28"
safe working load	600 lbs	600 lbs	600 lbs

NOA Twin Elite Riser

The NOA Twin Elite Riser™ Bed incorporates an extra-wide sleep surface for a new level of resident comfort and relaxation that makes them feel like they're back in their own bed. No more tossing and turning in standard-size long-term care beds.

WIDE SLEEP DECK

All NOA Twin Elite Riser beds come with an extra-wide 39" sleep surface.

Touch a Button

And watch the wide NOA Twin Elite Riser gently rise from 7.5" to 28" for resident safety and caretaker convenience.

Touch Another Button

And watch as the Twin Elite Riser adjusts the head and leg sections of the bed to fit your residents needs, and/or easily adjusts to a comfort chair position.

Guarantee a Good Nights Sleep

With 39" of sleep space your residents will have more freedom to move and turn in bed.



Features Include:

- A 39" wide sleep surface for superior comfort
- High-Low function from 7.5" to 28" for resident safety and caregiver convenience
- 600 pound safe working load
- Four quiet DC motors for smooth operation
- Effortless mobility in high or low positions
- Long linen life, excellent mattress aeration and easy cleaning
- Four Year Warranty on Electronics
- Lifetime warranty on bed deck and frame

Optional Features:

- Safety lock for immobilization
- Head and Foot Boards - Light Oak, Dark Oak, Wild Cherry, Walnut, Mahogany
- Soft Touch Side Rails
- Sub-acute handset
- Universal Assist Bar
- Floor Pad, 70" x 23.5" x 1"
- Therapeutic Foam Mattresses
- Footboard Mounted Staff Control
- Full Overbed Trapeze
- Deck Expansion System
- 42" Width Deck

ITEM	80" length	84" length
bed width	39.0"	39.0"
bed width, with side rails	44.0"	44.0"
overall bed length with head and foot boards	84.5"	88.5"
deck height range	7.5" to 28"	7.5" to 28"
safe working load	600 lbs.	600 lbs.

NOA Elite Bariatric Bed

Wide Sleep Deck

All NOA Elite Bariatric Beds come with a comfortable 42" wide sleep deck that easily adjusts from 80" to 84" in length.

Soft - Touch Side Rails

Side Rails are now made with a soft touch vinyl coated steel skeleton. The vinyl coating protects the residents from accidental bruising and the steel skeleton provides a rigid hand hold for ingress/ egress and repositioning.

Recliner Position

For added comfort, the bed can be put into a Recliner Position with the bed in a Reverse Trendelenberg position and the knee gatch and head deck raised for greater resident comfort.

Trolley System

NOA Elite Bariatric Beds come with a trolley system which allows movement of the bed from one room to the other or upright storage to save space when not in use.

NEW OPTIONS

- In Bed Scale
- Nurse Control Panel
- CPR Quick Release



Trolley System



Features Include:

- Soft-Touch Side Rail
- Polyethylene Head and Foot Boards
- Transport Trolley System
- T/RT Capable
- Baseboard Bumper
- Entrapment Compliant

Optional Features:

- Magnetic Lockout Handset
- Universal Assist Bar
- Floor Pad, 70" x 23.5" x 1"
- Therapeutic Foam Mattresses
- Footboard Mounted Staff Control
- Full Overbed Trapeze
- Deck Expansion System
- In Bed Scale with Digital Readout
- Underbed Lighting
- Drainage Hooks
- CPR Quick Release
- Nurse Panel

ITEM

bed width
 bed width, with side rails
 overall bed length with head and foot boards
 deck height range
 safe working load

80" length

42.0"
 47.0"
 84.5"
 9.5" to 30"
 800 lbs

84" length

42.0"
 47.0"
 88.5"
 9.5" to 30"
 800 lbs

NOA Elite R600

Self-Leveling

A built-in sensor in the control box communicates with the actuators to keep the bed level at all times.

600 lb. Safe Working Load

4 Year Warranty

on electronics and lifetime warranty on bed and deck frame.

Hospital Grade Casters

allow for easy transport within the facility.

Fully Adjustable

The Elite R600 provides trendelenberg positioning as well as cardiac chair positioning.



Features Include:

- Height Range : 9.5" at low position to 29.5" at high position
- Standard Hospital Grade 4" Locking Casters
- High density, white Polyethylene head and foot boards with handholds
- Rib Deck design for added durability and increased mattress ventilation
- Can be Moved in any position and at any height
- Goes Straight up and Straight Down

Options Available:

- Universal Assist Bar
- In-rail Pendant
- Lock Out Hand Set
- Floor Pad, 70" x 23.5" x 1"
- Therapeutic Foam Mattresses
- Footboard Mounted Staff Control
- Full Overbed Trapeze
- Standard Laminate Colors available
- Deck Expansion System
- Underbed Lighting
- Drainage Hooks
- Manual CPR Release Head Deck

ITEM	76" length	80" length	84" length
bed width	35"	35"	35"
bed width, with side rails	40.0"	40.0"	40.0"
bed width, with assist bar	36.5"	36.5"	36.5"
overall bed length with head and foot boards	80.5"	84.5"	88.5"
deck height range	9.5" to 29.5"	9.5" to 29.5"	9.5" to 29.5"


NOA Light

Able to effectively meet the needs of your clinical staff, the NOA Light provides a safe working load of 450 lbs. and offers a number of quality standard accessories.

NOA Medical's newest high/low bed system quickly rises from a low of 8.5" to a high of 25".

Fully adjustable

The four motors provides your residents with comfort positioning at the head and foot area, as well as height adjustment.

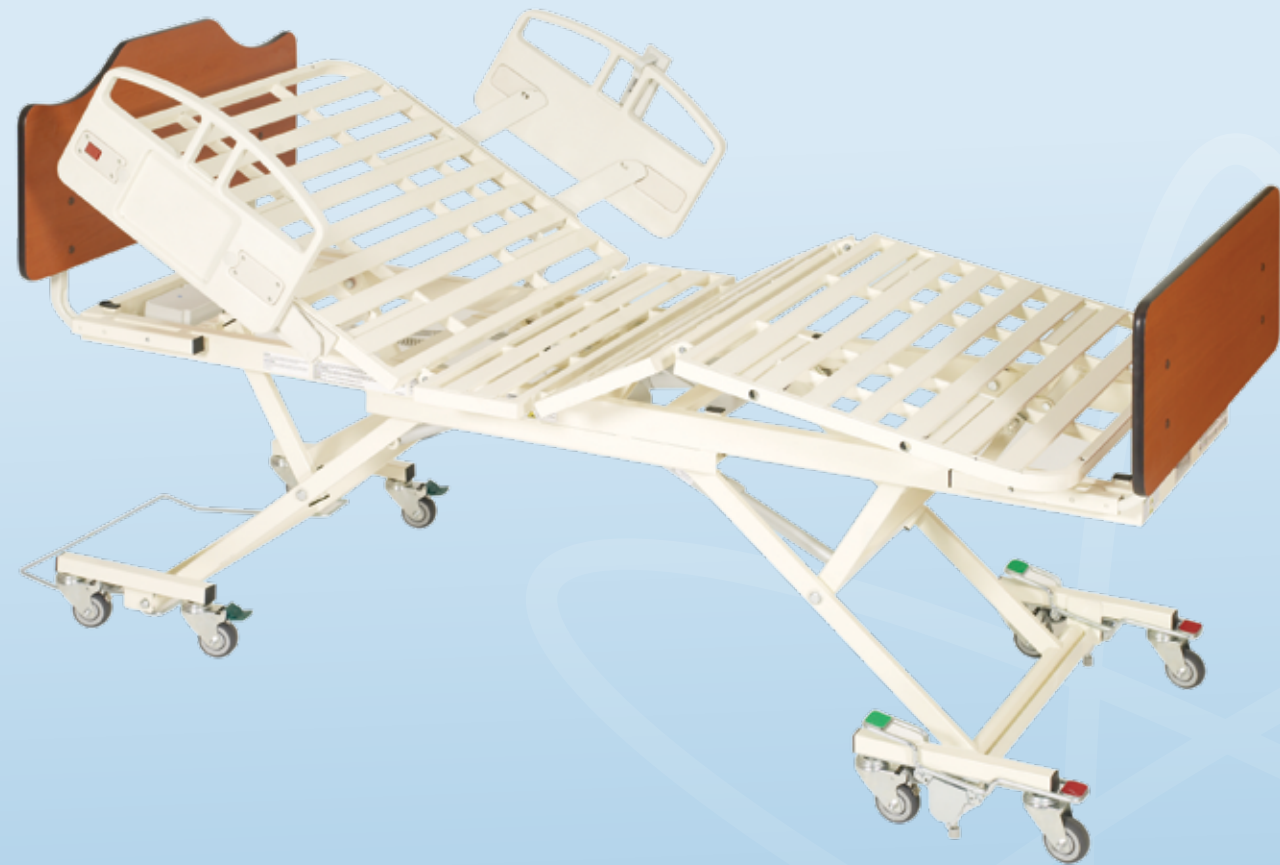

Features Include:

- Rises from a low of 8.5" to a high of 25".
- Has 4 DC Motors that provide elevation, head and foot adjustment.
- Has a sturdy rib deck providing a solid sleep surface and mattress ventilation.
- Is 35" wide by 80" long.
- Locking Casters.
- Can be moved in any position.
- Easy to clean and disinfect

Options Available:

- Soft touch side rails
- Assist bar
- Foot Board / Staff Controls
- Standard wood grain head and foot boards
- Deck Expansion System
- Safety Lock System

ITEM

80" length

bed width	35"
bed width, with side rails	40.0"
bed width, with assist bars	36.5"
sleep deck length	80.0"
overall bed length with head and foot boards	84.5"
deck height range	8.5" to 25"
safe working load	450 lbs

Accessories



Soft Touch Side Rails

Designed to provide extra safety for your residents while reducing the liability for the provider. Tall option also available.




Sub Acute Handset

Handset offers Trendelenburg/Reverse Trendelenburg positioning.




Standard Head and Foot Boards
(Available in 5 colors)

Durable melamine head and footboards are available in 36", 39" or 42" width.

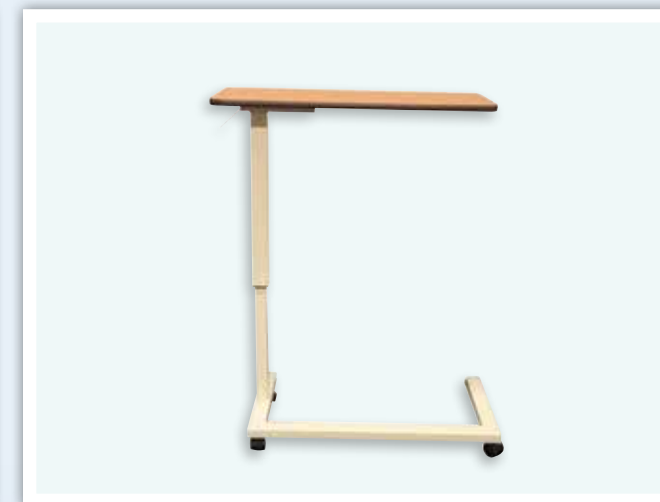

Over Bed Table
(Two models available: Gas Cylinder and Spring)

The NOA Overbed Table offers a solid, durable construction with a sturdy U-Base design. Featuring high pressure laminate tops available in our standard colors of Light Oak, Dark Oak, Walnut, Wild Cherry and Mahogany. Custom (special order) color choices are available.

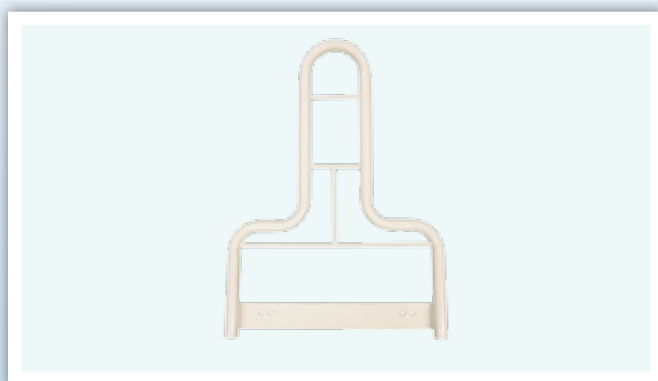

Standing Assist Bar



24" tall assist bar with a universal design mounts on either side of bed to provide residents with increased stability while repositioning, or entering and exiting the bed.



Folding Assist Bar

20" tall bar mounts securely to bed to provide increased stability for resident while repositioning, or entering and exiting the bed independently. This assist device easily swings out of residents way when not in use.

Color Options

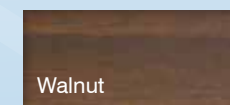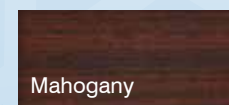

Accessories



Floor Mats

100% polyurethane mat designed with safety in mind. The beveled edge reduces trip hazards. Seamless design provides for easy cleaning and disinfecting.



Safety Lock

Provides immediate immobilization of the bed through an easy one step action. Easy to use.



In Rail Pendant

Eliminates searching for pendant control. Can be easily mounted on either side of bed rail and can face either resident or staff.



Scale

Bariatric model only

Available on the NOA Bariatric Bed. Easily weigh resident without having to move out of bed. Reads in Kg and Lbs.



Footboard Control


Provides the nursing staff with controls at the foot end of the bed. Provides easy lockout of all bed functions.



Under Bed Lighting

Available on the Elite Series of beds to provide a soft light to illuminate bedside and floor.

Accessories (cont.)

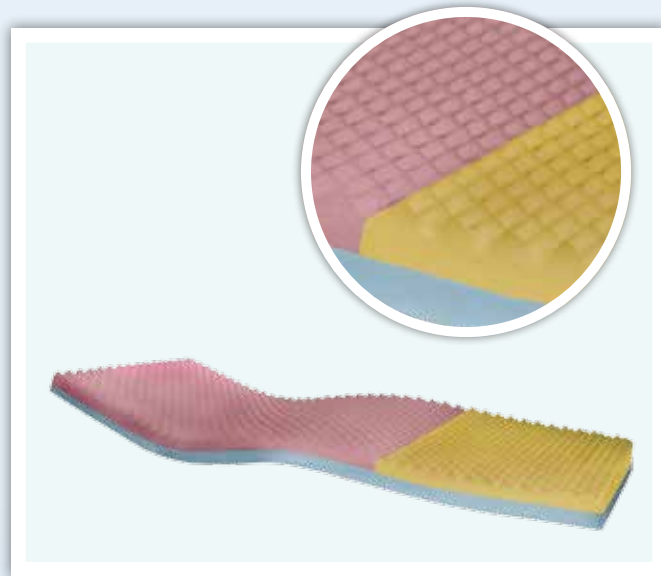

Mattress

Options available:

- Pressure Redistribution Foam Mattress
- Low Air Loss
- Alternating Pressure
- Custom Mattresses



Deck Expansion System

The Deck Expansion System allows you to grow a standard 35" deck to 39" and 42", our twin size decks from 39" to 42" and 45" and our bariatric bed system from 42" to 45" to 48". Our patent pending design allows for easy attachment to your NOA Bed Deck and offers you a quick and easy adjustment for your bed width requirements.



Trapeze

The NOA Trapeze is designed to help patients shift positions while in bed with minimal help.

- Provides safety and security
- Mounts directly onto bed
- Easy to assemble
- Economical
- Fits both standard and bariatric beds
- Rated at 800 lb

NOAH Hospital Bed Series



NOAH Hospital Platinum SC Bed

NOA Medical, through our new NOAH Hospital Division, introduces our newest product line for the acute care market, our NOA Hospital Bed Series. The NOAH Hospital Bed Series offers you three specific models to best fit your clinical and patient needs, while providing you with a combination of our standard and premium features, and price points to best meet your overall facility requirements.

Our top of the line NOAH Hospital Platinum SC Bed provides you with all of our standard features plus our premium bed package which includes a digital nurse panel with real time information and our in bed scale system. Our quality of product, progressive design and aggressive price makes our NOAH Hospital Bed Series one of the best values on the market.

Highlights

- Meets IEC 60601-2-52
- 600 lb. Safe Working Load
- Deck Height Range - 13" to 33.5"
- Self Leveling Motors
- 15 Year Warranty on Decks & Frame Welds
- 4 Year Warranty on Electronics
- Available in Both 80" and 84" Lengths



Electronic Features

- In Bed Scale with Digital Readout
- One Button CPR Quick Release Head Deck
- Out of Bed Alarm with Nurse Call Button
- Under Bed Lighting and Backlit Handset
- Battery Backup

Other Features

- High Impact Quick Release Head and Foot Boards
- High Impact Side Rails
- Central Locking System with 5" Casters
- Caster Base Covers
- Drainage Hooks

Controls

- Handset - 10 Button Backlit, Non-locking
 - » Head, Foot, Comfort Chair, Elevation, Trendelenburg, Under Bed Light
- Y-Cable (handset plugs into either side of bed)
- Nurse Panel
 - » All Handset Features Plus Angle and Height Indicators, CPR Button, Scale Readout
- Scale Control and Readout
- Pillow Speaker

Standard bed accessories



Under Bed Lighting

Night lighting to reduce falls and to assist in patient care.



CPR Quick Release

In the event of a CPR incident the head deck can be lowered in a controlled fashion by the push of the Electric CPR Quick Release Button on the Digital Nurse Panel. This feature provides quick patient access during an emergency.



Central & Directional Locking System

The central locking system locks the bed in place in one quick step. Central locking also provides for directional locking of casters for ease of transport.



5" Hospital Grade Casters

Hospital grade, anti-static caster provides for ease of movement and durability for long life.



Backlit Handset

Provides for ease of use for patients and staff at night.

Premium Bed Accessories



Digital Nurse Panel

Provides real time readouts of head and knee deck angle, T/RT angle, and bed height, and has a single push button for Quick Release of Head deck for CPR.

In Bed Scale with Digital Real Time Readouts

Weight displayed in either Lbs or Kg. Allows you to zero out for weight and add to while maintaining accurate patient weights.

NOAH Hospital Platinum NS Bed

Electronic Features

- Out of Bed Alarm with Nurse Call Button
- CPR Quick Release Head Deck
- Under Bed Lighting
- Battery Backup

Other Features

- Deck Height Range - 13" to 33.5"
- High Impact Quick Release Head and Foot Boards
- High Impact Side Rails
- Central Locking System with 5" Casters
- Caster Base Covers
- Drainage Hooks

Controls

- Handset - 10 Button Backlit, Non-locking
 - » Head, Foot, Comfort Chair, Elevation, Trendelenburg, Under Bed Light
- Y-Cable (handset plugs into either side of bed)
- Nurse Panel
 - » All Handset Features Plus Angle and Height Indicators, CPR Button
- Pillow Speaker



NOAH Hospital Platinum ML Bed

Electronic Features

- Out of Bed Alarm with Nurse Call Button
- Under Bed Lighting
- Battery Backup

Other Features

- Deck Height Range - 11" - 31.5"
- Manual CPR Quick Release Head Deck
- High Impact Quick Release Head and Foot Boards
- High Impact Side Rails
- Central Locking System with 5" Casters
- Caster Base Covers
- Drainage Hooks

Controls

- Handset - 10 Button Backlit, Non-locking
 - » Head, Foot, Comfort Chair, Elevation, Trendelenburg, Under Bed Light
- Y-Cable (handset plugs into either side of bed)
- Foot Board Control- all handset features plus battery power indicator
- Pillow Speaker



NOA

MEDICAL INDUSTRIES INC.

801 Terry Lane
Washington, MO 63090 **Toll**
Free: 1-800-633-6068
www.noamedical.com

About Us

NOA Medical Industries has been designing and manufacturing electric high/low beds and sleep safety accessories for over 30 years. Our progressive designs and attention to quality and function, combined with aggressive pricing, has made us one of the top providers of high/low beds to the long-term care marketplace. In 2014 NOA made a business decision to create a new division within NOA: NOAH Hospital Division. This new division introduced a line of hospital beds for the acute care market- the NOAH Platinum Series. This line includes dynamic features including an in-bed scale, central locking systems and a digital nurse control panel, to name a few. We approached it as a high functioning bed at a good price in order to help enter the market place.

As a small women owned business NOA is small enough to provide the personal attention you desire yet large enough to meet your needs in an appropriate time frame. Our well trained and experienced associates are empowered to ensure all NOA products meet our strict quality standards before leaving the factory. And, our warranty provides peace of mind for years to come.

NOA has 90,000 square feet of administrative, engineering and production space located in the heart of the country and serves customers across the U.S. and Canada.



FINANCING AVAILABLE

*Ask us about our
financing program*