



Platinum Behavioral Health Bed System

The NOA Platinum Behavioral Health Bed System provides a full electric, fully adjustable medical bed with the security features required for the safety of your patients and care givers.


Security Enhancements:

- All under bed and bed system wiring is secured and hidden when possible for added safety.
- Staff controls are embedded in foot end panel providing easy access for clinician. Staff control has lock outs to shut down individual settings/functions. A two-button code can also be used to shut down all functions while not in use. Light flashes when the security code is activated.
- Tamper resistant security screws and fasteners throughout the bed system.
- Bolted on mattress retainers.
- Central and Directional Locking with removable locking tool.
- HDPE secure deck covering for added safety for patient. Easy to clean/replace.



MEDICAL INDUSTRIES

Hoffmann Family of Companies

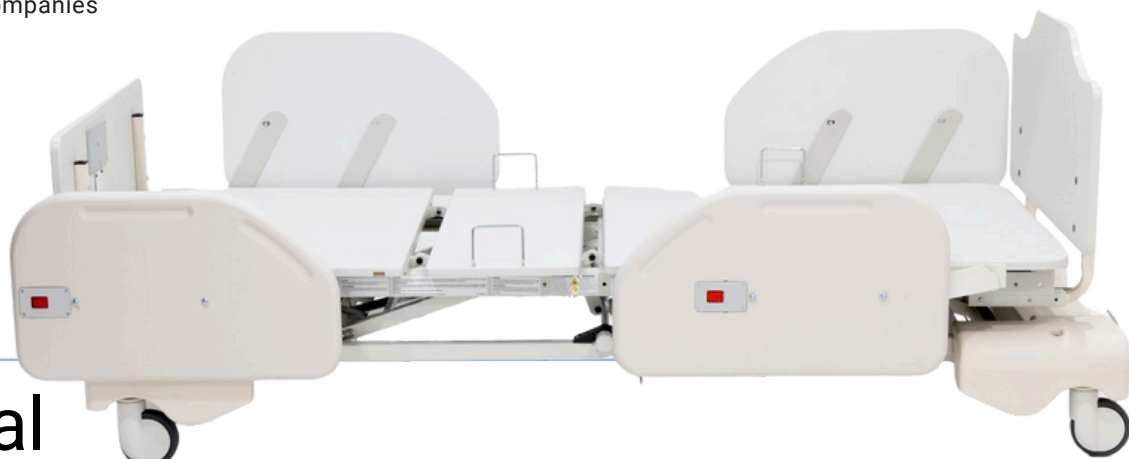

Technical Specifications

Platinum Behavioral Health Bed

Standard Features:

- HDPE, high impact head and footboards. Securely attached to the bed for added safety
- Heavy duty head and foot side rail – solid design with larger/deeper rail design. Securely fastened to the bed deck
- Removable power cord
- Central and directional locking system with removable locking tool
- 4 year electronic warranty and 15 year bed, deck, frame warranty

Features:

Dimension:

Overall Width (Includes Side Rails).....	39.5"
Overall Length (Includes Head & Foot Boards).....	86"
Patient Surface.....	34.5" x 80"
Bed Clearance.....	8"
Deck height range	17" - 32.5"

Optional Features:

- Under bed lighting
- Battery backup
- Zipperless mattress
- Removable patient hand pendent
- Patented integrated scale & bed exit alarm
- Lockbox

Features:

Dimension:

Safe working load	600 lbs
Trendelenburg & reverse Trendelenburg.....	14°
Max back angle	70°
Max hip-knee angle	35°

